### TECHNICAL SPECIFICATIONS



Handpiece: Not LED SP NEWTRON<sup>®</sup> (F12281) Irrigation: Connected to the water supply



Handpiece: LED NEWTRON<sup>®</sup> SLIM B.LED blue ring (F12900) white ring (F12905)

**Dental plaque disclosing liquid** F.L.A.G.™ for B.LED



Handpiece: LED NEWTRON<sup>®</sup> SLIM B.LED blue ring (F12900) white ring (F12905)

**Dental plaque disclosing liquid** F.L.A.G.™ for B.LED Irrigation: Connected to the water supply

Irrigation: 300 ml tank (500 ml tank in option: F62005)


Irrigation flow rate: 5 - 40 ml/min

NEWTRON P5XS BoleD® Handpiece: LED NEWTRON<sup>®</sup> SLIM B.LED blue ring (F12900) white ring (F12905)

Irrigation: 300 ml tank (500 ml tank in option: F62005)

Irrigation flow rate: 5 - 40 ml/min

**Dental plaque disclosing liquid** F.L.A.G.™ for B.LED **Bluetooth:** remotely sets power and irrigation

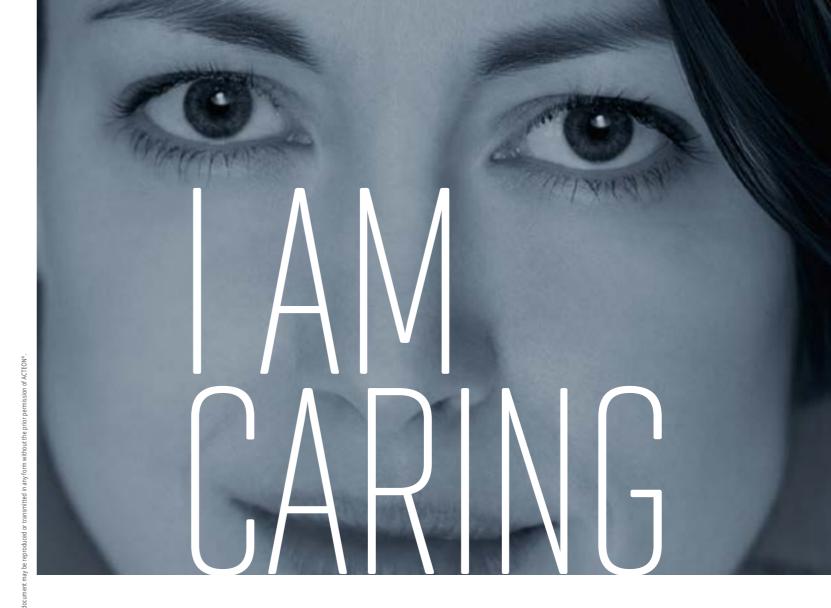



See individual product instructions manuals for further informations on **satelec.com/documents**.

SATELEC® S.A.S | A Company of ACTEON® Group 17 av. Gustave Eiffel | BP 30216 | 33708 MERIGNAC cedex | FRANCE Tel + 33 (0) 556 340 607 | Fax + 33 (0) 556 349 292 E-mail: info@acteongroup.com | www.acteongroup.com



### NENTRON

Preserve teeth with perfect ultrasonic vibrations and steel tip quality



## I TREAT ACCURATELY WITH PEACE OF MIND



### STEEL THE STANDARD

NEWTRON<sup>®</sup> tips hardness, designed with a specific steel, is the nearest to enamel.

With our tips, and their linear displacement, teeth are not only cared for but preserved.

3

# I CARRY OUT GENTLE TREATMENTS

#### **NEWTRON®** Technology, a guarantee of efficacy and safety

The NEWTRON® range ultrasonic generators are equipped with the exclusive and patented Newtron<sup>®</sup> technology, which brings to treatments more preservation, efficacy and comfort.

<b>IN</b> PRESERVATION	
<ul> <li>Controlled vibrations</li> <li>Steel tip quality</li> <li>Total irrigation control</li> </ul>	<ul> <li>More precise treatment and gentle on teeth</li> <li>Anatomy and enamel preservation</li> <li>No heating = safer for tissues</li> </ul>
<b>N</b> EFFICACY	
Frequency adjustment	<ul> <li>Maximal and continuous efficiency whatever the load applied</li> </ul>
Power regulation	→ Constant performances, no pressure needed
Powerful cavitation	<ul> <li>Fragmentation and elimination of deposits, cleaning and disinfecting effect</li> </ul>
<b>N</b> COMFORT	
Linear and regular vibrations	→ Tactile sense conserved
	→ Decreased fatigue for the clinician
	$\rightarrow$ Painless treatment for the patient




# I MEET ALL MY PATIENTS'

## CLINICAL NEEDS

#### Top of the range tips to provide first-class and caring treatments

- > The widest range in the market, having nearly 80 different tips, with exclusive designs, alloys and coatings for clinical versatility.
- ▶ World premiere: tips in **Pure Titanium for implant cleaning.**

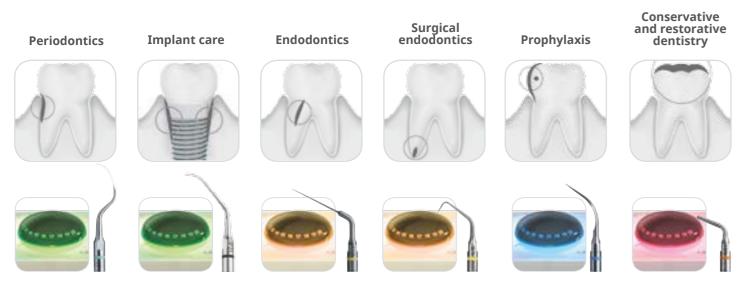

▶ Often copied, never equaled, our tips are conceived to vibrate with a single fundamental mode resonance, ensuring maximum and continuous efficiency.

#### **Intuitive settings for** secured clinical procedures

- ▶ Power use is accurately determined with the Color Coding System<sup>™</sup>, which intuitively associates each tip with one of the 4 power ranges, for safest interventions.
- ▶ The Bluetooth option on NEWTRON® P5XS automatically **configures power and irrigation** for a sustainable tip use, and maximum patient comfort. The **P5XS app** gives access to the tips range and proposes appropriate settings for each tip, as well as clinical videos and protocols.



### **Optimized external irrigation system**

- The 300ml or 500ml graduated tanks allow to mix solutions directly, and to fill in during procedures.
- ▶ The easy and precise flow adjustment provides more powerfull cavitation and maximum tip efficiency (deposits fragmentation, desinfecting effect).
- Any type of irrigation solution can be used including water, sodium hypochlorite or chlorhexidine.



## I SEE BETTER

## I TREAT BETTER

### A new and careful vision of scaling with B.LED technology

B.LED Technology guides practitioner in real time.

Applied on teeth, F.L.A.G.<sup>™</sup> for B.LED discloser targets dental plaque and turns fluorescent under the NEWTRON<sup>®</sup> SLIM B.LED handpiece blue light.

Exclusive benefits:

- ► Highlights dental plaque perfectly for a guided clinical care.
- Enhances treatment accurate healthy tissues.
- ► **Improves clinical quality** with more accurate and effective treatment.

This fluorescent revelation also e compliance.



Without B.LED

### VISUAL COMFORT QUALITATIVE TREATMENT ORAL HYGIENE EDUCATION

-WITRON BOLED




8

- ► Enhances treatment accuracy and avoids the overuse of instruments, thus preserving
- This fluorescent revelation also enables practitioners to educate patients and to encourage



With B.LED

# I CARE FOR DESIGN

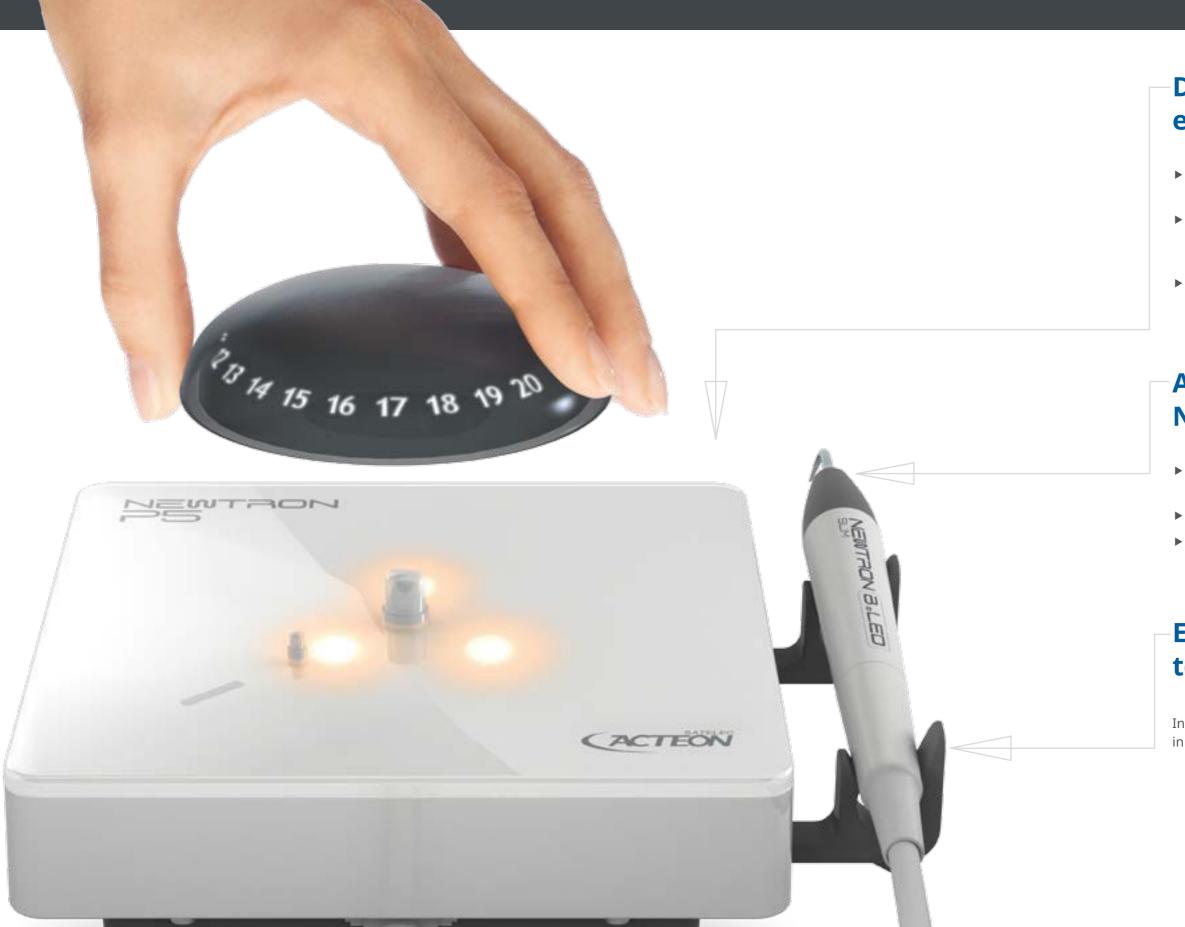

### Design responds to ergonomics and hygiene

- Elegant devices: flat glass surfaces, clean lines and luminous power dial.
- Adapted to practice: inclined front panel for better interaction with the practitioner and accessibility to the settings and the handpiece.
- Meeting hygiene requirements: total hygiene watertightness, removable power adjustment knob for easy decontamination.

#### A fine and balanced NEWTRON<sup>®</sup> SLIM handpiece

- Ergonomic shape regardless of the grip position, with flexible cord.
- ▶ Interior in Titanium for extreme resistance to corrosion.
- B.LED: dental plaque revelation and excellent intraoral visibility (2 white and blue interchangeable and autoclavable LED rings, and 6 High power LEDs).

### Exclusive handpiece holder to prevent falls

Installed on the front or the side of the device, the holders in silicone can be removed and autoclaved.